



# BG5 Function

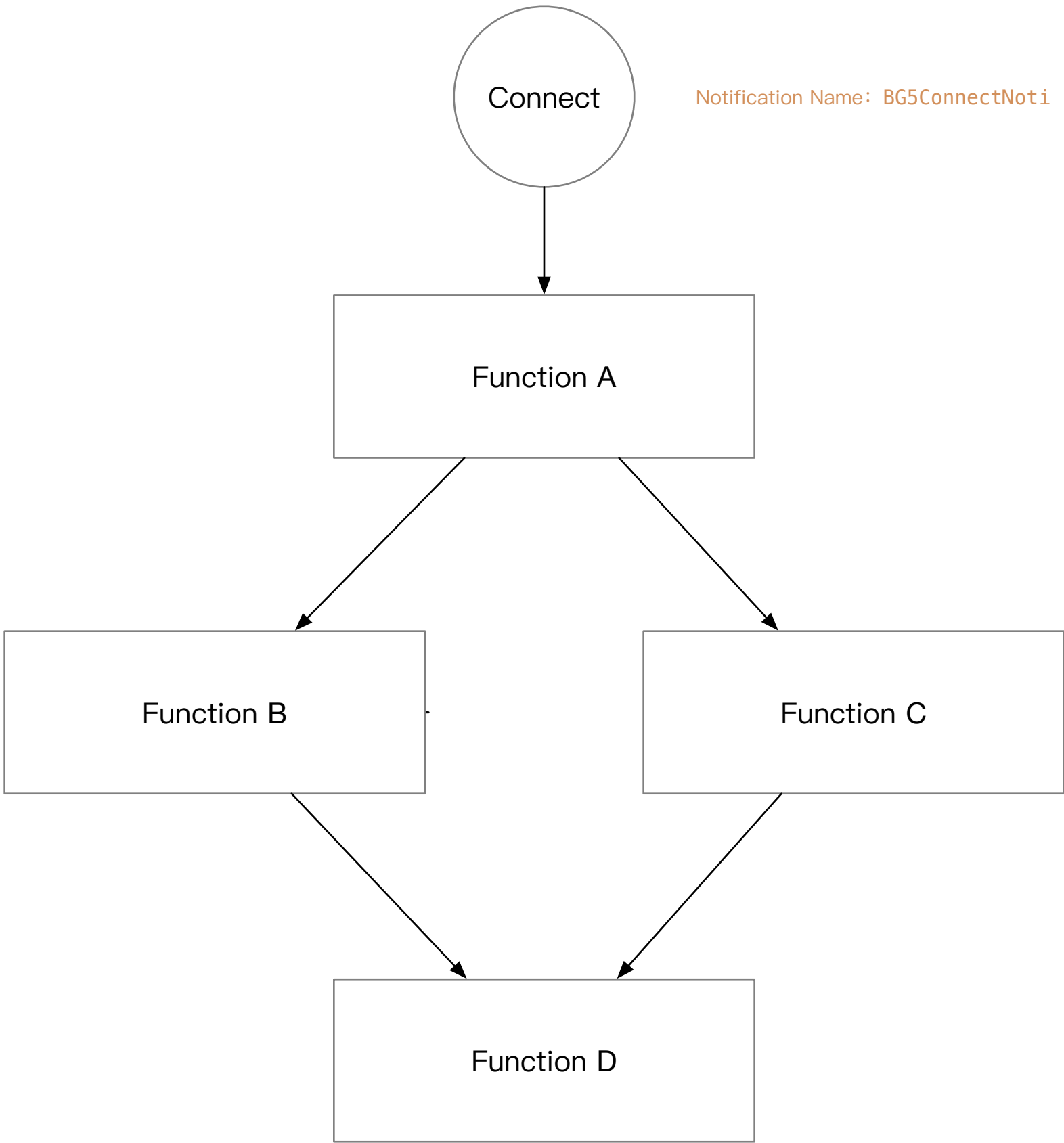

Note:

- 1、 API call is parallel call and sequential execution.
- 2、 Function B is for 1000I blood strips; Function C is for EGS2003 blood strips.
- 3、 Function D should be called after Function B or C.

# Function A: Set time, unit and sync offline data

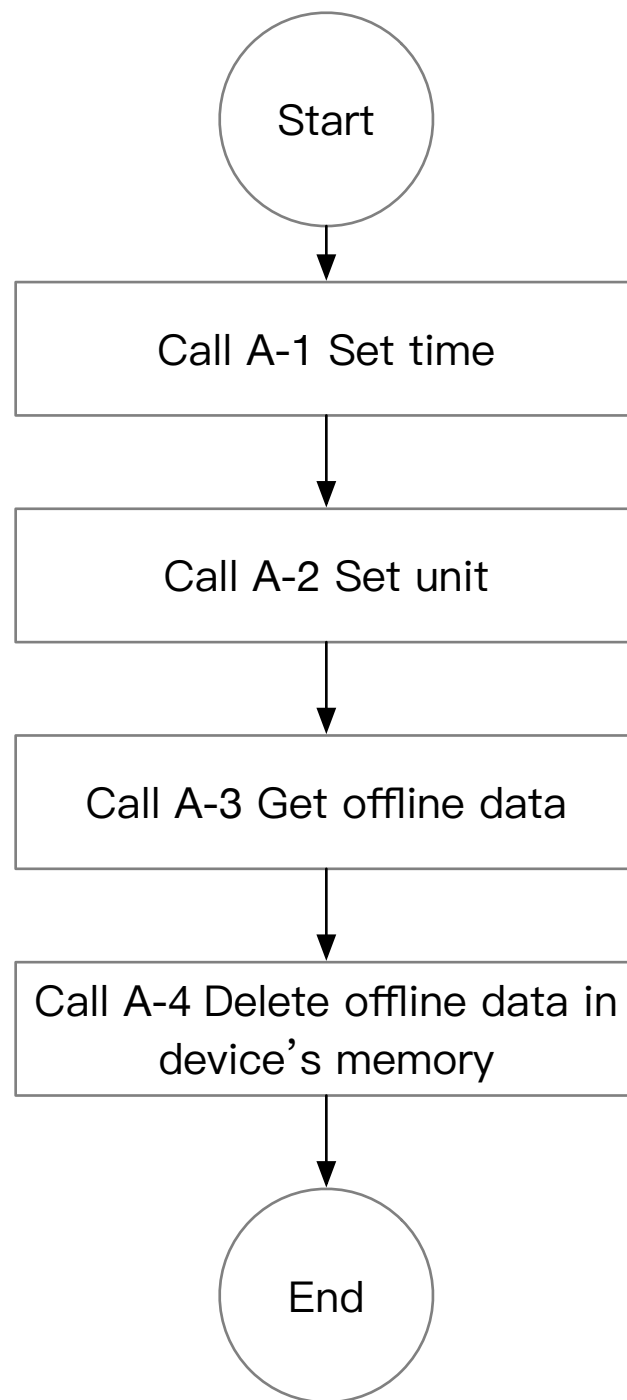

A-1

```
-(void)commandBGSetTime:(DisposeBGSetTime)disposeBGSetTime  
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```

A-2

```
-(void)commandBGSetUnit:(BGUnit )unitState  
    DisposeSetUnitResult:(DisposeBGSetUnit)disposeBGSetUnitResult  
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```

A-3

```
-(void)commandTransferMemoryData:(DisposeBGDataCount)disposeBGDataCount  
    DisposeBGHistoryData:(DisposeBGHistoryData)disposeBGHistoryData  
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```


A-4

```
-(void)commandDeleteMemoryData:(DisposeBGDeleteData)DisposeBGDeleteData  
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```

Note:

A-1 must be called before using the offline measurement, otherwise the time of offline data may be incorrect.

## Function B: Sending code flow for 1000I strips

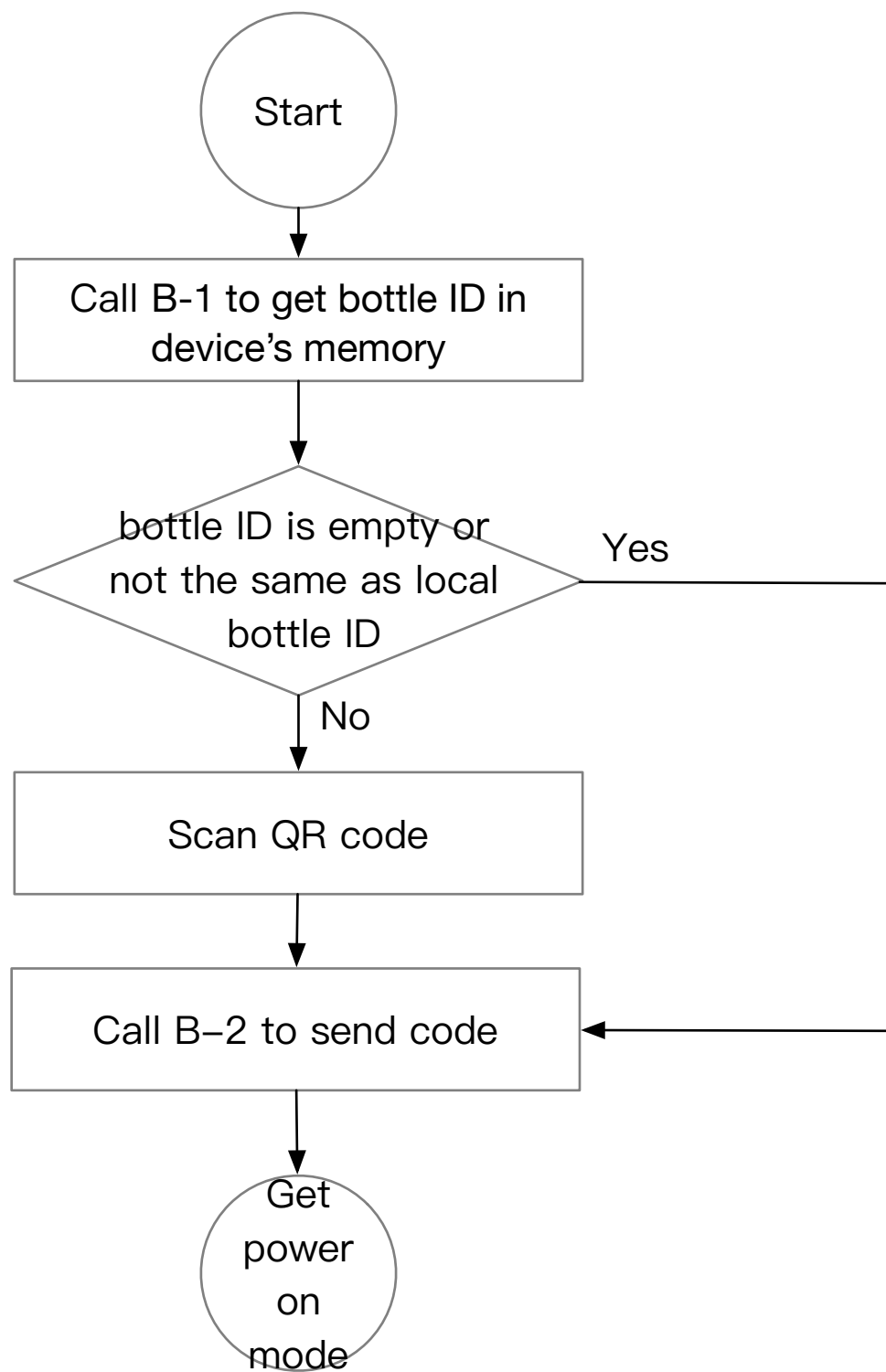


B-1

```
-(void)commandBGGetBottleID:(DisposeBGBottleID)disposeBGBottleID  
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```

B-2

```
-(void)commandSendBGCodeWithMeasureType:(BGMeasureMode)testType  
    CodeType:(BGCodeMode)codeType  
    CodeString:(NSString*)encodeString  
    validDate:(NSDate*)date  
    remainNum:(NSNumber*)num  
DisposeBGSendCodeBlock:(DisposeBGSendCodeBlock)disposeBGSendCodeBlock  
DisposeBGStartModel:(DisposeBGStartModel)disposeBGStartModel  
DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```


# Function C: Sending code flow for EGS2003 strips

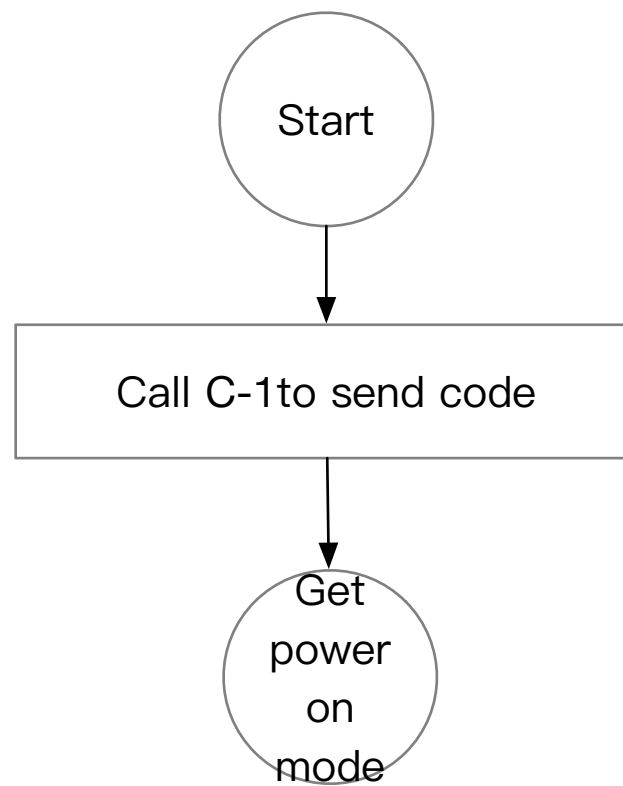

C-1

```
-(void)commandSendBGCodeWithMeasureType:(BGMeasureMode)testType  
    CodeType:(BGCodeMode)codeType  
    CodeString:(NSString*)encodeString  
    validDate:(NSDate*)date  
    remainNum:(NSNumber*)num  
DisposeBGSendCodeBlock:(DisposeBGSendCodeBlock)disposeBGSendCodeBlock  
DisposeBGStartModel:(DisposeBGStartModel)disposeBGStartModel  
DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```

# Function D: Online measurement flow

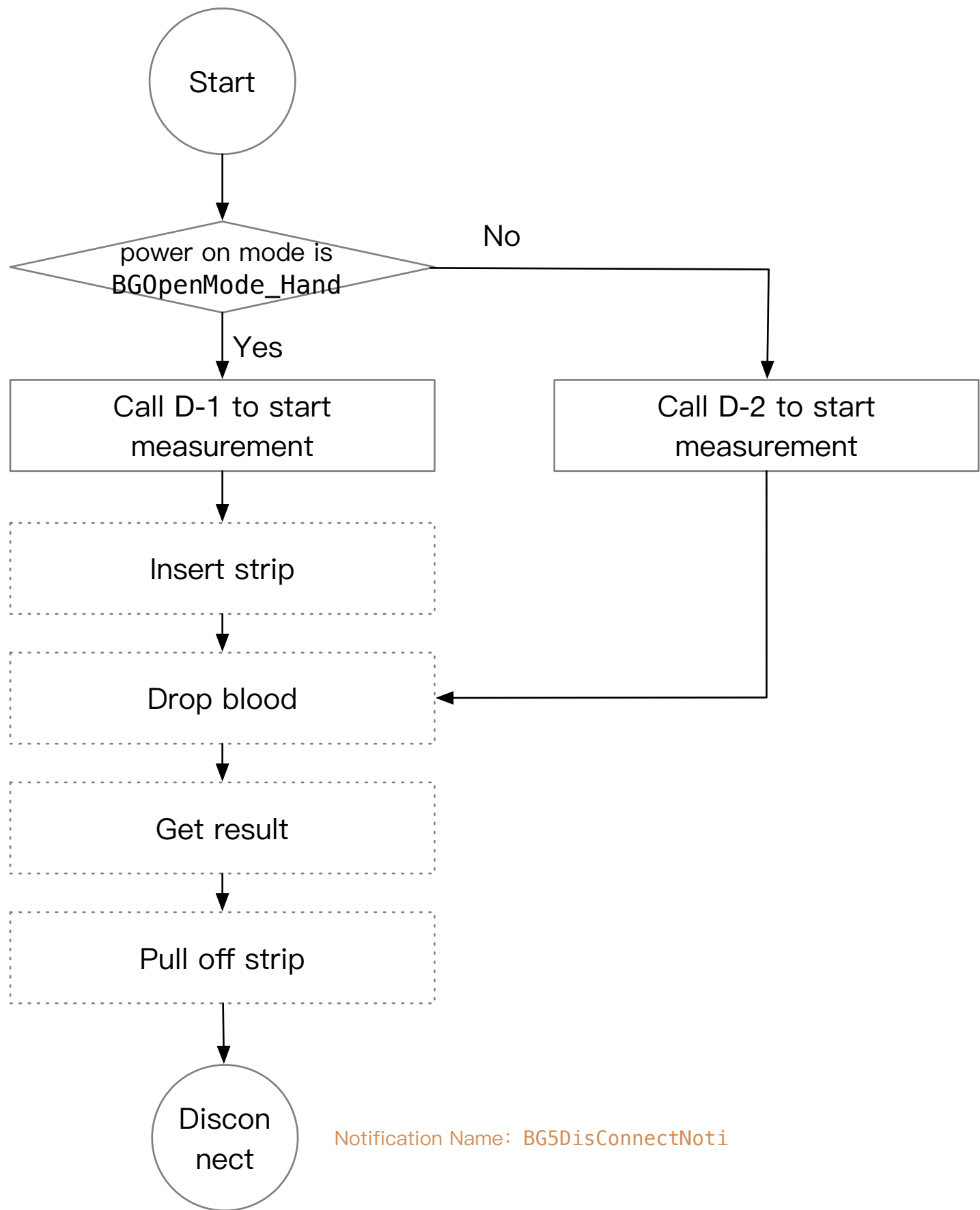

D-1

```
-(void)commandCreateBGtestModel:(BGMeasureMode)testMode
    DisposeBGStripInBlock:(DisposeBGStripInBlock)disposeBGStripInBlock
    DisposeBGBloodBlock:(DisposeBGBloodBlock)disposeBGBloodBlock
    DisposeBGResultBlock:(DisposeBGResultBlock)disposeBGResultBlock
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```

D-2

```
-(void)commandCreateBGtestStripInBlock:(DisposeBGStripInBlock)disposeBGStripInBlock
    DisposeBGBloodBlock:(DisposeBGBloodBlock)disposeBGBloodBlock
    DisposeBGResultBlock:(DisposeBGResultBlock)disposeBGResultBlock
    DisposeBGErrorBlock:(DisposeBGErrorBlock)disposeBGErrorBlock;
```